

# HemoCue® BT Connect



## Wireless Connectivity using Bluetooth® Technology

HemoCue® BT Connect offers wireless connectivity for HemoCue analyzers making it possible to send results from HemoCue analyzers to devices with Bluetooth 4.0 capability.

Enjoy true portability when needed. Connect to a PC or mobile device without using serial or USB cables. An excellent example of the use of HemoCue® BT Connect is the combination of HemoCue® Hb 201+ System, HemoCue® BT Connect and the Health-Trender Anemia solution for macro level trending of anemia prevalence.

With dedicated support and service together with training and education based on more than 35 years of experience, you can count on HemoCue for the right solutions for all your needs.

*Accuracy Starts With Us*

---

To learn more, please visit [hemocue.com](http://hemocue.com)

### Improve Your Workflow

- ▶ Send results wirelessly
- ▶ Portable connectivity

### Safeguard Patient Testing

- ▶ Avoid transcription errors with electronic transfer of results

### Easy to Connect

- ▶ Simple installation
- ▶ Multiple connector angles

### Ease of Use

- ▶ One button operated
- ▶ Clear LED indication
- ▶ Easy battery replacement

# HemoCue® BT Connect


## Product Specifications

<b>Weight</b>	Including battery: 26 g Excluding battery: 20 g
<b>Dimensions</b>	60 mm x 44 mm x 10 mm
<b>Storage Temp.</b>	-20 °C to +55 °C
<b>Operating Temp.</b>	+10 °C to +45 °C
<b>Relative humidity, non-condensing</b>	< 90%
<b>Battery Type</b>	CR2450, UL-rec.
<b>Wireless communication</b>	Bluetooth 4.1 -Bluetooth Low Energy Frequency Band: 2400 MHz – 2480 MHz Modulation Method: GFSK Max Output Power: 3dBm
<b>Interface</b>	Serial, 3.5 mm plug

## System Requirements

Device with Bluetooth® Wireless Technology 4.0/Low Energy (BLE) support  
 HemoCue® Hb 201+ Analyzer  
 HemoCue® Hb 301 Analyzer  
 HemoCue® Glucose 201+ Analyzer  
 HemoCue® Glucose 201 RT Analyzer

## Overview

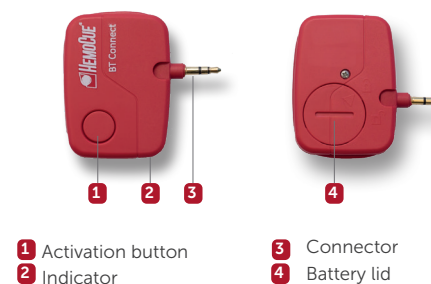

## Connect to the Analyzer



## Activate the HemoCue® BT Connect

Make sure the analyzer is on.  
 Press the activation button.



## Connector Angles

HemoCue® BT Connect has multiple connection alternatives:

Option 1



Option 2




Option 3



Option 4



## Indicator Status

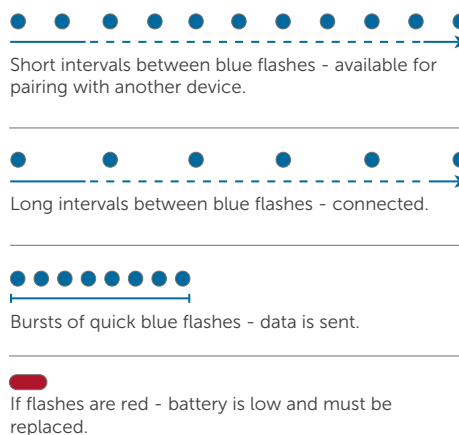

HemoCue AB | PO Box 1204 | SE-262 23 Ängelholm | Sweden  
 Phone: +46 77 570 02 10 | Fax: +46 77 570 02 12 | info@hemocue.se | hemocue.com

HemoCue has been a leader in point-of-care medical diagnostics for over 35 years. We specialize in giving healthcare providers lab-quality results comparable to that of a clinical lab.

